# Deluxe Adjustable Nasal CPAP Mask

Item# **CM106S** Small  
**CM106M** Medium  
**CM106L** Large

## Intended Use

The Sunset Deluxe Adjustable Nasal CPAP Mask is intended to be used in home, hospital, or institutional environments for patients who have been prescribed CPAP/VPAP therapy. This device is intended to be used under the specific direction of a physician.

Caution: Federal law restricts this device to sale by or on the order of a physician.

## Warnings

- This mask should only be used with CPAP systems recommended by your physician or respiratory therapist.
- Do not use this mask unless the CPAP or VPAP system is turned on and operating properly.
- Never block the mask's vent holes or exhalation port.
- Note: CPAP systems are intended to be used with special masks with ventholes in either the connectors or the mask that allows continuous flow out of the mask. When the CPAP machine is turned on and functions properly, new air from the CPAP machine flushes the exhaled air out through the mask's exhalation port. However, when the CPAP machine is not operating, enough fresh air will not be provided through the mask and exhaled air may be rebreathed. Rebreathing of exhaled air for longer than several minutes can in some circumstances lead to suffocation.
- If oxygen is used with the CPAP system, the oxygen flow must be turned off when the CPAP machine is not in operation. Note: When the CPAP device is not in operation and the oxygen flow is left on, oxygen delivered into the ventilator tubing may accumulate within the CPAP machines enclosure and create a risk of fire.
- Note: At a fixed flow rate of supplemental oxygen, the inhaled oxygen concentration will vary depending on the pressure setting, patient pattern and leak rate.
- Do not smoke if oxygen is in use.
- Contact a health care professional if chest discomfort, shortness of breath, stomach distention, belching, or severe headache are experienced.
- Caution: At low CPAP pressures the flow through the exhalation port may be inadequate to clear all exhaled gas from the tubing and some rebreathing may occur.
- Must only be used by persons who have received adequate training and instruction from qualified medical personnel.

## Contraindications


This mask may not be suitable for patients with the following conditions: impaired cardiac sphincter function, excessive reflux, impaired cough reflex and hiatal hernia. It should not be used if the patient is unresponsive, obtunded, unable to remove the mask, or uncooperative.

## Cleaning instructions


- Detach headgear from mask.
- Handwash gently with mild soap and warm water, rinse the mask and components thoroughly.
- Let mask and components air dry before reusing.
- Replace mask if there are cracks, tears or excessive discoloration.

⚠ Caution: The use of bleach, chlorine, alcohol or aromatic-based solutions may cause hardening of the materials of construction and reduce the life of the product.


## Instructions for use

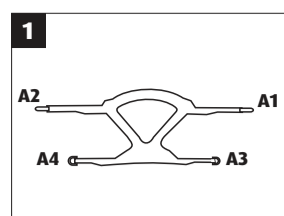
Place headgear on a flat surface with Velcro side down. Detach 2 Velcro strips on A1 & A2.

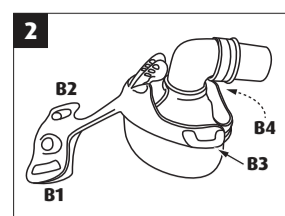
Identify headgear fastening positions on the mask.

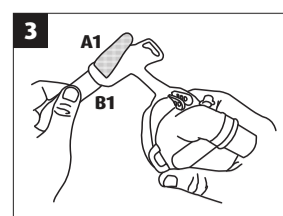
Pull the upper headgear strap A1 through the slot B1 on the mask, making sure the Velcro faces outward.

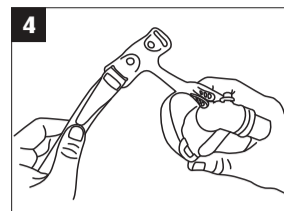Fasten the Velcro.

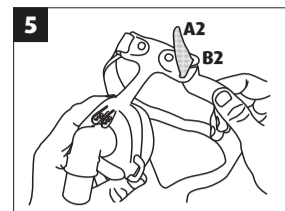
Pull the upper headgear strap A2 through the slot B2 on the mask, making sure the Velcro faces outward.

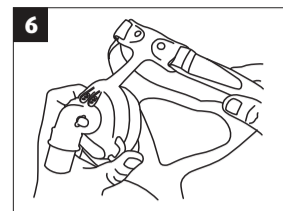Fasten the Velcro.

Place the headgear over your head and adjust both upper headgear straps to a comfortable tension.

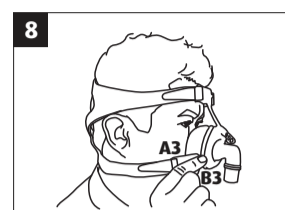
Hook the plastic buckle A3 onto B3 and A4 onto B4 on mask frame.

Readjust both upper headgear straps A1 & A2 to a comfortable tension. Squeeze forehead t-piece lock release C to tilt and adjust.

Adjust both lower headgear straps A3 & A4 to a comfortable tension.

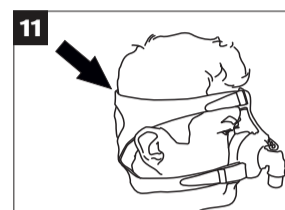Make sure the overhead strap is in proper position to avoid slipping down.

To remove the mask and headgear, simply unbuckle both buckles from the mask frame and pull the upper head strap away from the head.

To reuse the mask and headgear, simply place upper headstrap over head and hook up to both buckles.

Manufactured for  
 Sunset Healthcare Solutions  
 180 N Michigan Ave Ste 2000  
 Chicago, IL 60601

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