

Available from



# Immediate, simple treatment.



IT'S THE



By **AM** AIRWAY  
MANAGEMENT

Item#  
**CM039**

MAP: \$249

**Precision-Fit, Immediate,  
Simple and Affordable Snoring  
and Mild to Moderate Sleep  
Apnea Appliance.**

#### **WHY CHOOSE MYTAP:**

##### **IMMEDIATE TREATMENT**

Can be fitted in 15 minutes in a sleep lab, physician's office or HME.

##### **REMOULDABLE TRAYS**

Quick and easy to fit or refit. No bite registration needed.

##### **LOW PROFILE**

Small, low-profile, streamlined trays for patient comfort.

#### **WHO TO TREAT:**

Non-compliant CPAP patients  
Patients paying with cash  
Patients with high deductibles  
Patients not eligible for CPAP therapy



FLYMTP REV 1.01

**To order or for more information**

P 312.997.9980 E [sales@sunsethcs.com](mailto:sales@sunsethcs.com) W [SUNSETHCS.COM](http://SUNSETHCS.COM)





**The myTAP Concept**

The myTAP is a simple, patient-friendly treatment for snoring\* and sleep apnea. The design is based on the most effective oral appliance on the market, the TAP™, which is clinically proven in over 32 independent peer reviewed studies to treat snoring and obstructive sleep apnea (OSA). The unique trays of the myTAP are made from an innovative plastic that when heated becomes soft and can be molded to the teeth at the time of fitting. If reheated, the trays will return to their original shape and can be refitted.

**How myTAP Works**

The myTAP advances and stabilizes the mandible preventing the soft tissues of the throat and tongue from collapsing into the airway, which causes snoring and sleep apnea. The adjustment mechanism is attached to the lower tray and is threaded through the adjustment post of the upper tray. An adjustment dial is screwed onto the end of the threaded adjustment mechanism allowing the patient to adjust the mandible forward. The myTAP hardware allows the patient to control the treatment position with incremental advancements as needed.

**What is "Precision-Fit?"**

"Precision-Fit" is a category for the myTAP trays created by Airway Management. The accurate molding capabilities of the trays coupled with the unique ability for the trays to be remolded if necessary makes the myTAP superior to non-custom oral appliances. It also features thin, sleek trays which create a comfortable, low profile appliance.


**Why Choose the myTAP?**

- Precision-Fit trays can be remolded or refitted as necessary to achieve the perfect fit.
- The adjustable design creates active patient participation in therapy, promoting compliance and success.
- Midline advancement provides the patient the ability to achieve the ideal treatment position.
- The simple design and ease of use makes it an affordable option to treat snoring,
- An AM Aligner comes with every myTAP to help offset potential bite changes.
- The myTAP can be titrated in the sleep lab.
- Patient can titrate the myTAP at home and then return for a follow up sleep study.

**How does myTAP help grow my business?**

By offering the myTAP, you can tap into a pre-existing group of patients who need an effective treatment solution.

- Treatment of snoring and sleep apnea is one of the fastest growing markets.
- 1 in 4 of the population\* has some degree of sleep disordered breathing.
- Approximately 50% of diagnosed patients are non-compliant or forgo traditional treatments (like CPAP) altogether.
- Can be an alternative treatment for CPAP users when traveling.
- Affordable treatment for patients paying with cash or who have a high deductible.
- A way to treat patients who changed providers due to Medicare contracts.

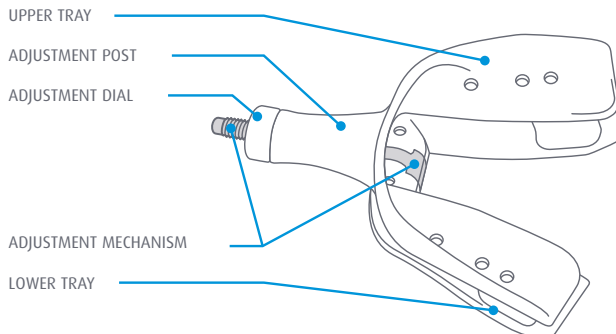

Snoring may be an indicator of a more serious medical condition called Obstructive Sleep Apnea (OSA). The patient should seek medical advice if OSA is suspected.  
 \*Based on the US adult population

**To order or for more information**

P 312.997.9980 E sales@sunsethcs.com W SUNSETHCS.COM

