

BONGO^{Rx}

A Sleep Apnea Therapy Device

EPAP AIRAVANT
MEDICAL

Item#

CAPBNG500


Bongo Rx all sizes starter kit

- All four sizes: small, medium, large, and extra large
- Travel case
- Drying stand
- Headgear

- No machine, hose, electricity, or batteries needed
- Use at home, traveling, camping... *use anytime, anywhere*
- Small enough to fit in a shirt pocket
- Reusable and easy to clean with simply soap and water
- Designed, molded, and assembled in the USA
- Clinically tested to demonstrate safety and effectiveness
- Supplemental lifestyle option for your compliant patients
- Patient friendly option for your non-compliant patients

The Bongo Rx is a simple device with soft nasal seals and is available in various sizes (S, M, L, and XL). The Bongo Rx is an EPAP (expiratory positive airway pressure) device and works as follows:

- During inhalation, small valves open to allow normal breathing through the device.
- During exhalation, the valves close and direct exhaled air through small vents to generate EPAP, which keeps the airway open.

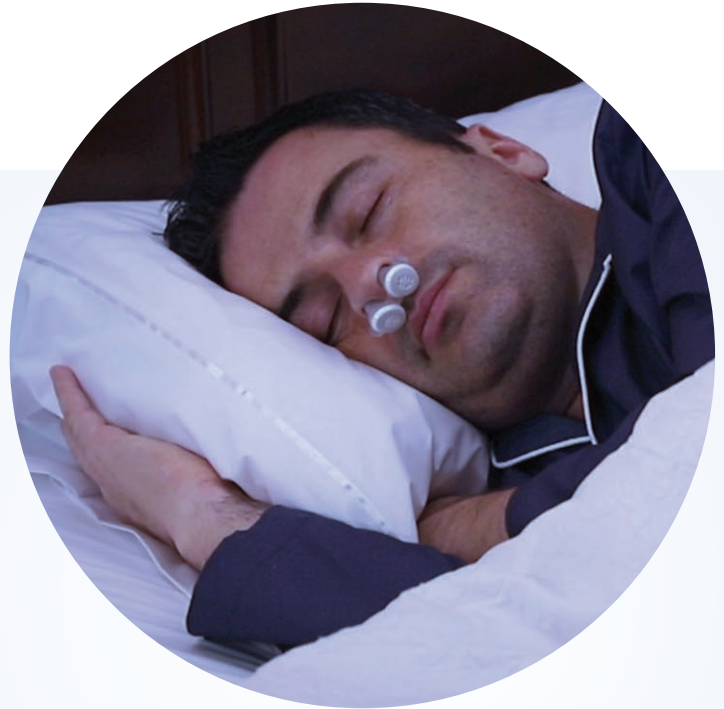

Annual Replenishment Pack

Includes 4 pairs, travel case, drying stand, and headgear

- CAPBNG502** Small
- CAPBNG503** Medium
- CAPBNG504** Large
- CAPBNG505** Extra Large



Single Replenishment Unit

Includes 1 pair and headgear

- CAPBNG512** Small
- CAPBNG513** Medium
- CAPBNG514** Large
- CAPBNG515** Extra Large